

## Series parameters

	3L Standard	5L Standard
Capacity	3L	5L
Product size	55*36.5*28cm	55*36.5*28cm
Chamber size	175*55*280 (L*H*Dmm)	220*68*320 (L*H*Dmm)
Packing size	64*44*38cm	64*44*38cm
GW/NW	40/33KG	40/33KG
Vacuum Pump	Diaphragm Type	Diaphragm Type
LCD screen	Keybutton control	Keybutton control
Single tray	(28.6*18.5cm) with caradle	(31.6*20.5cm) with caradle
remarks		



For Details Contact :

**MEDISTER**  
Infection Control Systems Pvt. Ltd.

Regd. Off. : G-119, Shagun Arcade, Premises, Film City Road, Gen. A. K. Vaidya Marg, Near Dindoshi Depot, Malad (East), Mumbai - 400 097.

Factory : Unit R6-R6A, Qureshi Indl. Area, Sagar Industrial Estate, Zilu Parab Rd, Chinchoti, Off Kaman Bhiwandi Naka, Naigaon (E), Tal.:Vasai  
Phone: +91-8291040850 / +91-8291029242  
E-mail : admin@medister.in • Website : www.medister.in

Our Sales / Service Network :


Maharashtra / New Delhi / West Bengal / Tamilnadu / Karnataka / Andhra Pradesh / Kerala / Madhya Pradesh / Gujarat / Rajasthan / Punjab / Uttar Pradesh / Orissa / Assam / Himachal / Uttarakhand / Goa / Bihar / Jharkhand / Nagpur / Chattisgarh.



**Mini&Flash**  
Autoclave

# 20 cycles records

When printer without paper, the system will keep 20 cycles records.

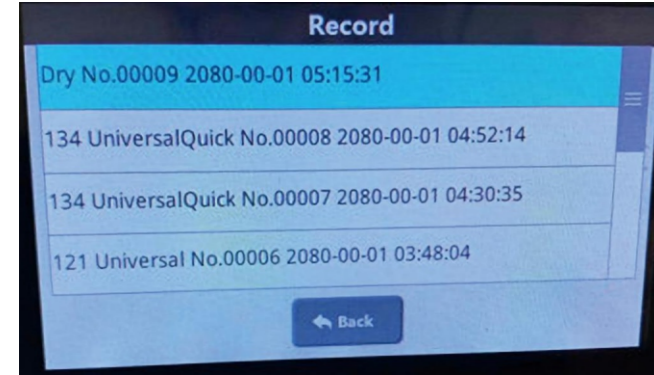

# Touch Screen Operation System

Easy & Friendly





# Fast Sterilization Time

134°C Solid 10mins (8mins sterilization + 2mins dry)



# Better Dryness

Aluminum alloy Rectangle Chamber  
With better dryness than stainless steel Cylindrical Chamber  
Individually Dry Cycle (5mins)

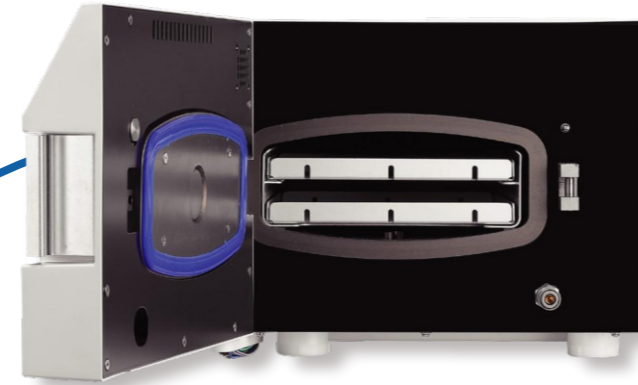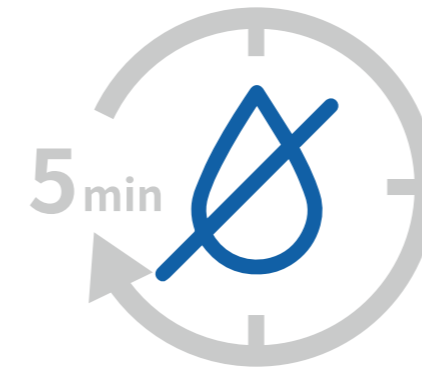

# More secure



**1 Chamber Test**  
The chamber we use which also reliable after 100 thousand pressure test(6KG).



**2 Sealing Ring**  
Special design of sealing ring which can stand the test of 3KG pressure.



**3 Double security door lock**  
Latch with diameter 12cm; locks hook with two grooves.



**4 Short circuit protective device**  
Shut down the machine in case of electricity leakage or short circuit.



**5 Safety valve**  
Prevent the chamber from over pressure.